

Endoscopy

CleanFlow™ Insufflation Filter

Bacterial and Viral Protection

The hazards of surgical smoke have been recognized by AORN, OSHA, NIOSH and ANSI. The CleanFlow insufflation filter prevents bacterial and viral contamination from possible backflow of surgical smoke, protecting both equipment and operating room personnel.

Advantages

- Virus Filtration Efficiency (VFE) of 99.999% for viruses to 0.02 microns.**
- Hydrophobic membrane provides complete fluid backflow protection.
- Meets all ECRI and OTT recommendations.

Maximum Flow

- Compatible with insufflators with flow capacities to 24 liters/minute.

* Data on file

+ Bacterial and viral retention at the flow rates up to 24 liters/minute.



medovations®

Specifications:

Flow Rate: 24 liters/minute
Bacterial Retention: 99.999%
Virus Retention: 99.999% to 0.02 microns*+
Hydrophobic: Hydrophobic to 1 psi
Latex Free, Sterile

Model	Description	Count Per Case
SU1101	Insufflation Tubing with CleanFlow Filter	10
SU1102	Insufflation CleanFlow Filter	10
SU1114	Insufflation Tubing with CleanFlow Filter and CPC Connector	10
SU1122	Insufflation Tubing with CleanFlow Filter—22mm Connector	10



SU1101



SU1102



SU1114



SU1122

Endoscopy • Liposuction • Neurosurgery



medovations, inc. 102 east keefe avenue, milwaukee, wi 53212 usa
800.558.6408 414.265.7620 fax 414.265.7628 medo@medovations.com www.medovations.com

VersaProbe™ Suction/Irrigation and Veress Needle



The Versatile Solution


The patented VersaProbe adapts to most irrigation systems including bags, bottles, pumps and cuffs and allows universal integration of various probes and/or electrodes. Its soft rear port connector permits utilization of electrodes and laser guides.

Reliable Irrigation

VersaProbe's precision-engineered trumpet-valve design provides high suction, clog resistant aspiration and leak-free irrigation. A comfortable saddle grip design enables left- or right-hand use, and our custom engineered tubing is memory-free for better drapability.

Cost Effective

VersaProbe converts into an electrocautery suction/irrigator to compete against higher priced systems.

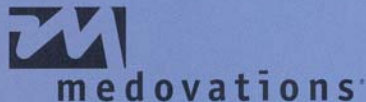

Model	Description	Count Per Case
SU2001-2S	5mm/33cm Probe with 2-spike Irrigation Tubing, sterile	12
SU2001-2P	5mm/33cm Probe with Dual Suction Pump Irrigation Tubing, sterile	12

Hi-Flow Insufflation Needle

Ergonomically designed non-slip handle. Blunt tip, spring-activated stylet with visual abdominal free-space indicator. Precision, sharp needle for ease of use, 120mm and 150mm lengths, disposable.

Model	Description	Count Per Case
SU1390	120mm Insufflation Needle, sterile	20
SU1395	150mm Insufflation Needle, sterile	20



Endoscopy • Liposuction • Neurosurgery